

iSnare[®] system




listening... and delivering solutions[®]

Maintain scope positioning...

The iSnare® system from US Endoscopy is a family of unique, multi-functional devices ideal for endoscopic mucosal resection and saline assisted polypectomy.

The iSnare® system's unique design...

- Features alternating injection needle and polypectomy snare operated by a multi-functional handle
- Facilitates immediate injections in response to bleb absorption/dissipation or bleeding
 - Minimizes wasted procedural time – inject or snare at any time throughout the procedure without device exchange while maintaining scope position
 - Features a spring-loaded, luer lock handle for consistent needle projection and retraction

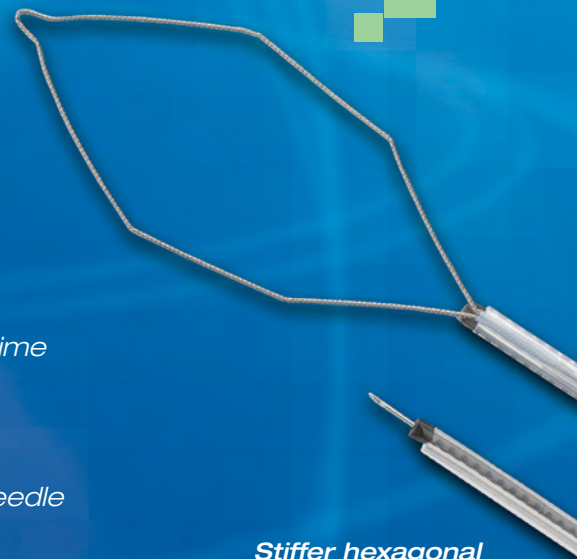

Stiffer hexagonal snare style

The iSnare® system can be utilized for all of the following endoscopic needs:

SAP (saline assisted polypectomy), EMR (endoscopic mucosal resection), hemostatic injection and adherent clot removal, injection therapy, tattooing pre- and post-polypectomy, cleansing, applying dye or stain during chromoendoscopy, and any other procedure requiring injection therapy and snare resection techniques.

For more information on the iSnare® system, contact US Endoscopy today.

product number	description	needle gauge	needle length	snare size	sheath diameter	sheath length	active cord connector	units/box
00711085	oval snare	23	5mm	2.5 x 4.0cm	3.0mm	230cm	Olympus	5
00711086	oval snare	25	5mm	2.5 x 4.0cm	3.0mm	230cm	Olympus	5
00711088	hexagonal snare	25	5mm	2.5 x 4.0cm	3.0mm	230cm	Olympus	5
00711089	hexagonal snare	23	5mm	2.5 x 4.0cm	3.0mm	230cm	Olympus	5

Visit www.usendoscopy.com to see product videos, literature, and much more.



listening...and delivering solutions®

5976 Heisley Road
Mentor OH 44060

phone 440 / 639.4494

fax 440 / 639.4495

customer service 800 / 769.8226

www.usendoscopy.com

760205 Rev. F